


## ABOUT HIGH FOCUS SHOT PLUS

### PRODUCT INFORMATION

#### Operation Instructions

- Press and hold the ⑦ Power Button for 2-3 seconds to turn on the device. When the power is on, the ⑧ LCD Screen will light up with a notification sound and the ② LED Indication will automatically turn red.
- After the device is turned on, it will be automatically set to Level 1 with a voice prompt. Press the ⑦ Level Button once to adjust the level. Each time you press the ⑦ Level Button, the level will increase by 1 and the newly set level will be displayed on the ⑧ LCD Screen. \*There are 5 intensity levels. When using the product for the first time, be sure to start from Level 1.
- Once the level is set, ensure that the ③ Skin-boosting Head and ④ Ultrasonic Head are in full contact with the skin.
- Press the ① SHOT Button to begin delivering automatic shots. \*Please use with the dedicated gel applied, and continue moving the device to avoid delivering continuous shots to the same spot. \*A notification sound is provided for each shot. \*A notification sound is provided every 100 shots.
- After 600 shots, an additional care standby screen will be displayed on the ⑧ LCD Screen for 10 seconds. If there is no button/app operation within the 10 seconds, the device will automatically shut down. If you wish to continue the additional care, press the ① SHOT Button once within 10 seconds. \*To turn off the device during operation, press and hold the ⑦ Power Button for more than 2 seconds.
- For clean and safe storage, always clean the ③ Skin-boosting Head and ④ Ultrasonic Head with a dry cloth or cotton pad after use. \* Be careful not to damage the black film on the head.
- To charge the device, connect the included Type-C charger to the ⑤ Charging Port located at the bottom rear side of the device. \*During charging, the charging status will be displayed on the ⑧ LCD Screen along with a charging connection sound, and the device will not operate while charging.

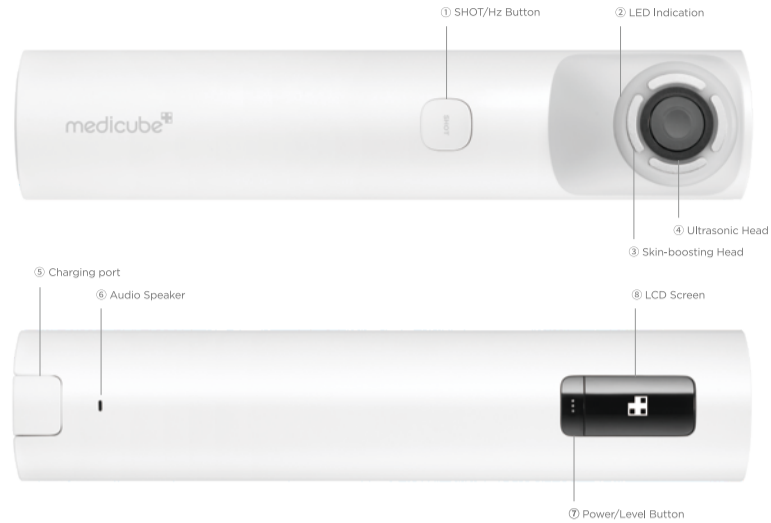

04

### CHECKLIST FOR BLUETOOTH

- Remote Control  
LED Color, volume changes and skin-boosting on/off can only be made within the remote control screen after Bluetooth has been connected. The volume level and LED color set with the remote control are saved as default values, and when the device is restarted, the LED color for the next session will automatically be set to the last color used.
- Automatic Reports  
Usage data recorded while the device is not connected to the app is temporarily stored in the device and automatically uploaded to the app's report once connected. The usage levels and LED color are saved in the report based on the longest-used level and color.
- Multiple Bluetooth Devices  
Only one device can be connected via Bluetooth at a time. If you wish to connect to another device after Bluetooth pairing, please turn off the power of the currently connected device before attempting to connect to another device.
- Pairing Failure  
Please ensure that the device is not already connected to another smartphone. If the device is connected to another smartphone, open the app on the connected smartphone, tap on the [HIGH FOCUS SHOT PLUS Banner] or [Connect Bluetooth Banner], then navigate to [Settings] → [Disconnect Device] to disconnect the connection. Afterward, try reconnecting again.
- Switching between Smartphones  
Open the app on the connected smartphone, tap on the [HIGH FOCUS SHOT PLUS Banner] or [Connect Bluetooth Banner] then navigate to [Settings] → [Disconnect Device] to disconnect the connector. Afterward, try connecting to the new smartphone.



### PRECAUTIONS

#### Precautions During Use

- Before use**
- Do not irritate the skin before use
  - This device is intended for use on the face and neck. Do not use it on other body parts.
  - Stop using the device if the black film on the head has been damaged.
  - If the product is not operating correctly, do not disassemble it; instead, contact customer support.
  - Do not tap the head or the black film that contacts the skin directly on hard surfaces. This may cause malfunction or damage.
  - Remove all metal accessories in contact with the user's body.
  - Do not use the device for purposes other than its original intended use; it may lead to accidents.
  - Do not use the device in a bathtub, shower, or water-filled reservoir
- During Use**
- Apply a sufficient amount of the dedicated gel to the desired area.
  - Ensure that the head makes contact with the skin.
  - Do not deliver shots if the head is not in full contact with your skin or if it's not touching your skin at all.
  - Avoid using the device on the following areas: between the eyebrows, around the eyes, eyeballs, philtrum, lips, around the mouth, nose, jawbone, Adam's apple, thyroid area, wounds or areas where procedures or surgeries have been performed.
  - \*The manufacturer and seller are not responsible for any side effects caused by using the device on prohibited areas.
  - If you experience excessive heat or discomfort, immediately remove the device from the skin.
  - Move the device continuously over the skin during operation without pausing over any one area.
  - When using the product for the first time, always start at the lowest intensity level and adjust as needed.
  - Follow the recommended usage frequency (once or twice a week). Otherwise, it may

- cause skin irritation.
- The device will not operate while charging.
- The appliance is not to be used by individuals (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction.
- Children should be supervised to ensure that they do not play with the appliance.

#### After Use

- You may experience skin irritation after use depending on your skin condition. If you feel any irritation, discontinue usage.
- Avoid any activity that may irritate your skin after use. This includes using hot tubs or saunas.
- Avoid using any exfoliating products for 7 days after use.

#### Storage/Handling Precautions

- Please follow the usage and storage conditions below.

Operating Temperature	5°C to 35°C
Operating Humidity	0% to 90%
Storage Temperature	1°C to 60°C * Please note that storing the device at temperatures below 0°C may cause freezing.
Storage Humidity	0% to 90%

- Always store the device safely after use to prevent damage to the black film.
- Store in a place out of the reach of children and away from moisture.
- Do not store in places where magnetic fields or electromagnetic waves are generated; it may cause malfunctions.
- After use, clean the head with a wet wipe, clean gauze or soft cloth and store it dry.
- Ensure that the device's power is turned off before storage.
- Do not use flammable substances such as paint thinner or benzene for cleaning, as it can cause chemical reactions, discoloration, and deformation.
- Be careful not to drop or shock the device, as it may cause malfunction or damage.
- Ensure that no water or cosmetics enter the inside of the main unit.
- Do not use the device in conjunction with other electrical appliances such as electric heating pads.

05

- Do not use the device in hot places such as near fire or direct sunlight, and do not place the main unit or cable near heating appliances.
- Be cautious when using the device as it may malfunction if it comes into contact with other devices.
- Do not strike or drop the main unit on hard objects.
- Avoid using the device in high temperature, high humidity, or dusty environments, such as bathrooms, kitchens, and hot cars.
- When using the device, maintain a distance from electric appliances such as other beauty devices, TVs, and radios.
- Do not use the device at voltages other than the rated voltage. Using it where the rating is exceeded can generate heat, possibly leading to fires, electric shock, or malfunction.
- Do not handle the device or attach/detach the charging cable with wet hands; it can lead to electric shock.
- Do not attempt to repair, disassemble, or modify the main unit yourself; it can cause malfunctions.

#### User Precautions

If you are in poor health, under the influence of alcohol, taking sleeping pills, feeling fatigue, or under medical treatment, please do not use the device. If you experience pain or unusual symptoms during use, discontinue use immediately. **Repeated use of the device on the same area or improper usage may cause damage to the skin. If any abnormality appears on the skin, cool the affected area and consult a medical professional immediately.**

\*Even in normal use, individuals with sensitive skin may experience slight redness. Please read the instructions and precautions carefully and perform proper aftercare.

#### Do not use the device if you have the following conditions:

- Individuals with symptoms of redness and stinging due to excessive sun tanning
- Individuals with allergies, especially those with sensitive skin
- Individuals using electronic devices worn directly inside the body, such as cardiac pacemakers
- Individuals using electronic medical devices for life support, such as artificial heart-lung machines

- Individuals using implanted electronic medical devices such as electrocardiogram devices
- Individuals using steroid-based hormone medications for an extended period
- Individuals with liver function disorders who have symptoms of vascular dilatation syndromes in the circulatory system
- Individuals undergoing outpatient treatment
- Areas with inflamed acne
- Areas with dermatitis caused by cosmetics, or individuals with atopic dermatitis
- Areas with blemishes accompanied by itching, or areas with blemishes caused by physical irritation
- Areas with metal, plastic, silicone, or other foreign objects inside the body
- Patients with blood disorders, heart diseases, acute illnesses, bleeding areas, malignant tumors, fractures, infectious diseases, eyelids, trauma areas, skull areas, and other children, pregnant women, breastfeeding women, and elderly individuals
- Individuals having paralysis or difficulty communicating
- \*This product is intended for cosmetic use only. It is not intended to diagnose, treat, cure or prevent any medical condition or disease.

### SPECIFICATION

**Model Name** ME-HFSP **Product Name** MEDICUBE AGE-R HIGH FOCUS SHOT PLUS  
**Seller** APR Co., Ltd.  
**Seller Address** 36th floor, 300, Olympic-ro, Songpa-gu, Seoul, Republic of Korea  
**Customer Service** +82-1577-0719  
**Manufacturer** APR FACTORY Co., Ltd  
**Manufacturer Address** 50, Wonseokjeongyongso-gil, Poseung-eup, Pyeongtaek-si, Gyeonggi-do, Republic of Korea / 107, Gasan digital 2-ro, Geumcheon-gu, Seoul, Republic of Korea

**KC Registration Number** R-R-pgl-ME-HFSP  
**Instruction and Precaution** Refer to user manual  
**Serial Number and Manufacturing Date** Refer to manufacturing label

**Product Specifications**  
 Product Size 38.0 X 186.0 mm  
 Input Voltage : DC 5V / max. 2A (USB-C)  
 Battery : 3.7Vdc / 2,550mAh (Li-ion Battery)  
 Power Consumption : max. 30W  
 Weight : 280g  
 Charging Time : Approximately 3.5 hours  
 MADE IN KOREA

**Battery Specifications**  
 Quality Indication in accordance with the Electrical Appliance and Consumer Products Safety Control Act Safety Confirmation Report No. ZU10601-25001

**Model No.** ND18650-2600mAh-3C **Product Name** Li-Ion Battery  
 Rated 3.7Vdc / 2,550mAh **Manufacturer** DONGGUAN LEAN POWER NEW ENERGY TECHNOLOGY CO.,LTD. Made in China

06

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:  
 (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:  
 - Reorient or relocate the receiving antenna.  
 - Increase the separation between the equipment and receiver.  
 - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected. Consult the dealer or an experienced radio/TV technician for help.

FCC Note Caution: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

### EXCHANGE / RETURN POLICY

#### Warranty

For all Age-R devices, we offer a 1-year limited warranty period from the date of your original purchase. If your device does not work properly due to a manufacturer's defect, we will replace the defective device free of charge within 1 year. This warranty does not cover the devices damaged by:

- Accident, misuse, or alteration
  - Use with unauthorized accessories or use other than instructed
  - Connecting to incorrect current and voltage (We only provide Type-C USB cable)
- This warranty only extends if the devices are used in conjunction with authorized device accessories (such as the enclosed charging cables).

#### Returns

For the Age-R devices, all returns must be initiated within 30 days of receiving your order. Please refer to the information as shown below. **Simple Remorse:** A refund due to a change of mind is only applicable within 30 days of receiving and will only be accepted if the product is returned in original, sealed, and unopened condition. \*In the case of a return for simple remorse (change of mind), the customer needs to bear the shipping fee. For international orders: You will be responsible for paying for customs fee and returning shipping fee. Additionally, we cannot provide free return shipping labels for international customers. **Damaged or defective product:** 1) Please reach out to us via our Customer Service Center email with the following details: First name and last name, email address, phone number, order number. Kindly attach a photo or video describing your issue. 2) We will proceed to issue a confirmation email for you (please note that our Customer Service usually takes 2-3 business days to respond). 3) After receiving a confirmation email, please send the product back to our warehouse via the return address. You will be informed of the return address after the email confirmation. 4) Once the product has arrived at our warehouse, our staff will inspect the condition of the product. 5) In case of product defect, we will proceed to exchange and send a new product to you.

\*Note: This process may take up to 1 month, including the time taken to ship the product back to our warehouse, and the time taken to inspect the product.

#### Subject of Recommendation of Lithium Battery Discarding and Charging Method

Warning: Since it may accompany danger such as heating, fire, and explosion, please observe the following matters. **1) Handling Matters** Since a swollen battery that can be visually identified may be dangerous, immediately contact the manufacturer or the seller. Use only the designated charger. Keep away from fire (Do not put inside a microwave oven). Do not leave it inside a car in summertime. Do not store or use at high temperature and humid place such as sauna. Do not use it for a long time while placing it on blanket, electric pad, and carpet. Do not store at a sealed space while leaving the power on. Be cautious not to contact the electrical terminal with metal products such as necklaces, coins, keys, and watches. Use the secondary lithium battery warranted by the manufacturer. Do not implement behaviors such as disassembling, compressing, and penetrating. Do not apply abnormal impact such as dropping it from a high place. Do not expose to a temperature above 50 degrees Celsius. Do not contact with humidity. **2) Disposal Guideline** Do not arbitrarily dispose used up battery. When disposing the battery, do not incinerate or put it inside a fire. Explosion and fire may occur due to heat. The method of battery disposal may vary according to countries and regions. Dispose with appropriate method. **3) Recommendations regarding Charging Method** Do not charge the battery with the method other than the one shown in the user manual. Fire or explosion may occur.

**Complete Charge Time : About 3.5 hours** \*Charge time may vary depending on the battery remnant amount of portable battery or surrounding temperature or charging environment. It may change without notice for performance improvement in specifications and exterior appearance.

themedicube.co.kr

Copyright © 2025, APR Corp.  
 All Rights Reserved.

07

